


Versatile


"EGG" slit

Digital Camera Compatible L-0229-E


L-0229-E Digital Camera Compatible Slit Lamp Microscope "EGG" is ideal for taking digital images both still and movie by our digital video camera beam splitter L-0541DC-1.3. All the extended functions were thoroughly designed for your comfort.

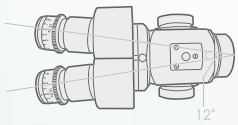


(Light source side)

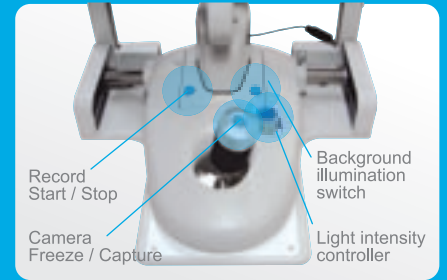

The Joy-stick features a shutter release, a video recording function, a background illumination switch and a light intensity control. Everything is on the unique shaped base of this slit lamp for ease of use. Still images and videos can be taken in one position.

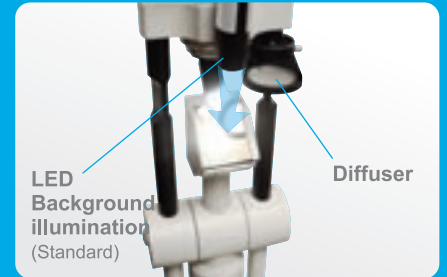

L-0229-E features the LED background illumination system which allows you to take clear and bright images of slit light. The light intensity of LED is set within the extent that does not give too much damage to human eyes. The light diffuser is useful when the full view of an eye is observed.

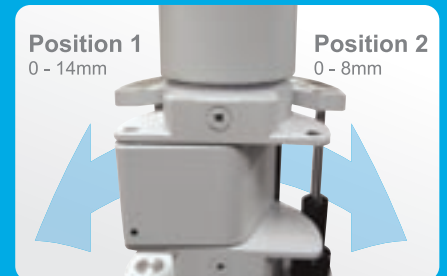

The unique function of L-0229-E called "S.I. GRANDE" enables you to optically change the illumination field up to 14mm by single click action. This is ideal for contact lens fitting observation. "S.I. GRANDE" is registered as Inami's patent.

SPECIFICATIONS

Microscope Head Type	Convergent Type Binocular Microscope
Magnification Changer	5 Steps by Drum rotation
Stereo Angle	12°
Magnification Ratio (Field of View)	6x (38mm) 10x (24mm) 16x (15mm) 25x (9mm) 40x (6mm)
Slit Width	0-14mm
Filters	Cobalt blue, Heat Absorption, Red-Free, 1/2 ND, Yellow
Weight	23 kgs
Dimensions	550(W)x400(D)x690(H)mm