


L-5110 / L-5130

APPLANATION TONOMETER



INAMI Applanation Tonometer L-5110 / L-5130

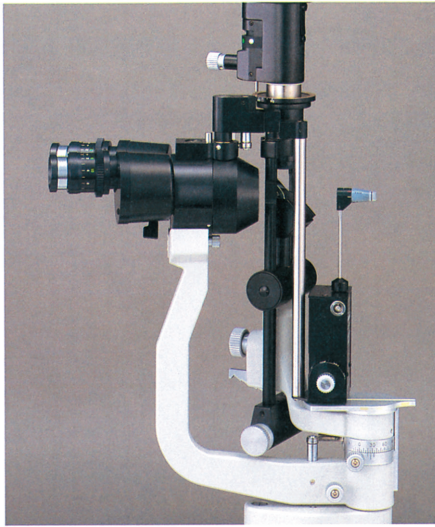

L-5110 (Swing arm type) mounted on INAMI Slit lamp

Swing arm Type

L-5110 INAMI Applanation Tonometer for the measurement of intraocular pressure is designed to be mounted on the microscope permanently.

When the tonometer is in use it can be swung forward just in front of the microscope and exact centering of the prism with the left eyepiece is engaged by notch. Illumination through the built-in cobalt filter projected to the plastic measuring prism for flattening the vertex of cornea gives reflection-free observation of split images.



L-5130 (Plate type) mounted on INAMI Slit lamp

Plate Type

L-5130 INAMI Applanation Tonometer can also be used with various INAMI Slit Lamps.

L-5130 can be placed on the guide plate in one of two possible positions. These positions are connected with the microscope, for observation can be made either through the right or the left eyepiece.



L-5110 / L-5130 contains 2 cone prisms.

Specifications:

Measurement range : From 0 to 80mmHg in 2mmHg increments $\pm 0.5\text{mmHg}$

Diameter of the pneumatic face : 3.06mm

Measurement : 47mm wide x 30mm deep x 85mm high

Weight : 725gs.

Standard Set

Applanation tonometer x 1

Prism x 2

Calibration bar x 1

Adaptor x 1

Specifications are subject to change without prior notice.

A Passion for a Nice View



Inami & Co., Ltd.

No.24-2, Hongo 3-chome, Bunkyo-ku,
Tokyo 113-0033, Japan
Tel: +81 (3) 3814-5916 Fax: +81 (3) 5684-2126
E-mail: corpinfo@inami.co.jp

