

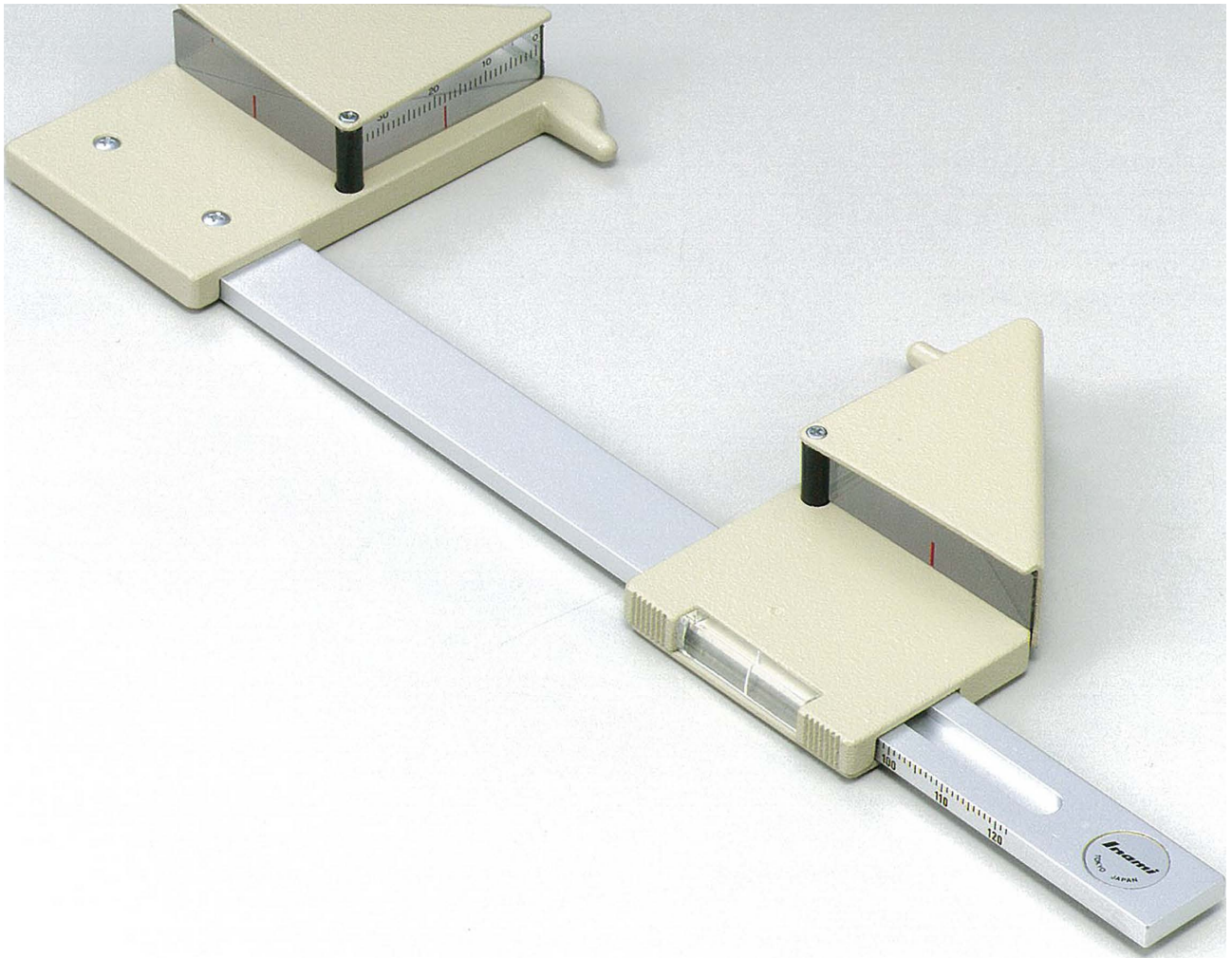


Inami
EXOPHTHALMOMETER
K-0161

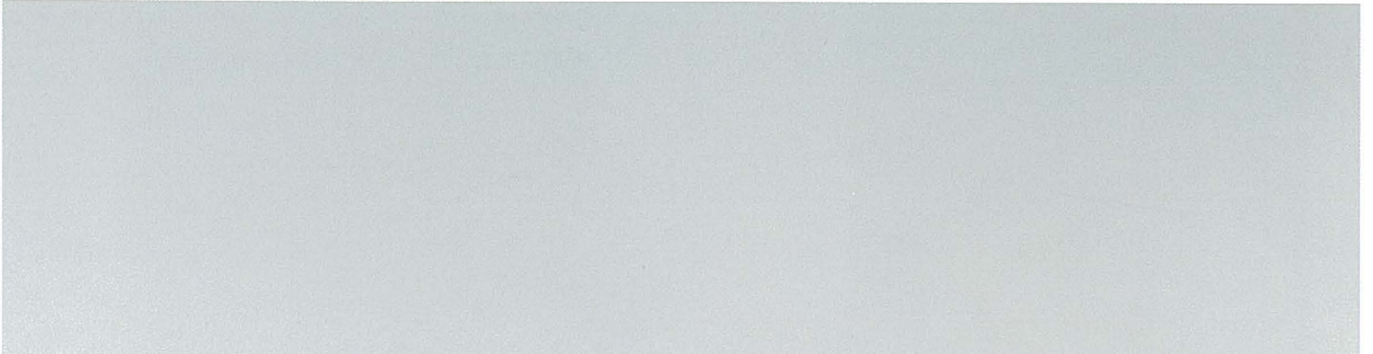

INAMI EXOPHTHALMOMETER, HERTEL K-0161

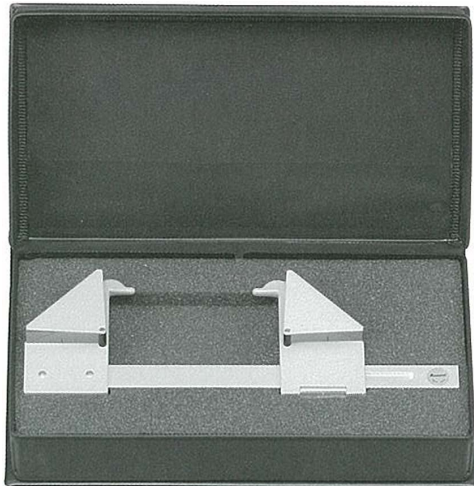

The Inami Exophthalmometer (Hertel) allows for precise measurement of corneal projection (free from parallax error) through simple alignment of twin prism, linear vertical targets. When these targets are overlapped, the degree of projection can be read directly off the gauge (mm).

■ Specifications

Weight:	350g
Length:	230mm
Case:	112(W)×265(L)×32(H)mm

● Specifications subject to change without prior notice.



Inami & Co., Ltd.

No.24-2, Hongo 3-chome, Bunkyo-ku,
Tokyo 113-0033, Japan
Tel: +81 (3) 3814-5916 Fax: +81 (3) 5684-2126
E-mail: corpinfo@inami.co.jp

EU REP
WMDE B.V.
Bergeweg 18
6085 AT Horn
The Netherlands

